TRANSFER PATIENTS WITHOUT Manual Lifting





Step1

Place the chair by the bed in the stretcher position.



Step 2

Attach the bed sheet to the transfer bar.




Step 3

Transfer the patient by turning the handle.



Step 4

Incline the chair to transport patient.

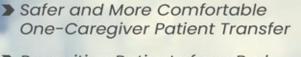

Reposition Patients from Bed to Chair Position WITHOUT Manual Lifting

Home Health and Hospice Approved



Patient Requirements

Fully

Accredited

Medicare

Provider

- The Barton Chair multi-positional patient transfer system is covered by Medicare when transfers between bed and a chair, wheelchair, or commode are required, because without the use of a lift, the beneficiary would be bed confined.
- Require supine/lateral transfers
- Non-weight bearing and chairbound patients that require Max-Assist transfers



Qualifying Patient Diagnoses

- Alzheimer's/Dementia
- Amputations (AKAs and BKAs)
- Multiple Sclerosis (MS)
- Joint Contractures
- Decubitus Ulcers (Stage II or Greater on Trunk or Pelvic Region)
- Degenerative Joint Disease (DJD)
- Osteoarthritis
- Paralysis (Paraplegia and Quadriplegia)
- Rheumatoid Arthritis
- Seizure Disorder
- Stroke

Medicare Covered and Coded E1035 & E1036



BARTON CHAIR Patient Transfer System

H-250

Weight Capacity: 250 lbs.



Transfer Patients WITHOUT Manual Lifting

I-400

Weight Capacity: 400 lbs.



Extended Width and Weight Capacity to Comfortably Accommodate Larger Patients

I-700

Weight Capacity: 700 lbs.



Extended Width and Weight Capacity to Comfortably Accommodate Larger Patients

Convertible Chairs	Part no.	Seat width	Overall width	Overall length	Overall height	Maximum user weight
H-250 Chair	B1002-BL	17"/ 43 cm	24"/ 61 cm	Prone: 72.5-80"/ 184-203 cm Wheelchair: 49-67"/ 124-170 cm	Prone: 27"/ 68 cm Wheelchair: 53-60"/ 135-152 cm	250 lbs / 113 kg
I-400 Chair	B1041-TS-GE	L 21"/ 53 cm	30"/ 76 cm	Prone: 72.5-78.5"/ 184-199 cm Wheelchair: 45.5-62.5"/ 115-159 cm	Prone: 27.5"/ 70 cm Wheelchair: 51-56"/ 130-142 cm	400 lbs / 181 kg
I-700 Chair	B1071-TS-GE	L 25.5"/ 65 cm	35"/ 89 cm	Prone: 72.5-78.5"/ 180-199 cm Wheelchair: 45.5-62.5"/ 115-159 cm	Prone: 27.5"/ 70 cm Wheelchair: 51-56"/ 130-142 cm	700 lbs / 318 kg
Patient Transfer System Part no.		Included in kit				
Patient transfer system kit for H-250		b2004e	2 chair brackets, 5 transfer straps, 2 bed posts, 1 mattress cover (non-slip), 1 Transfer bar, 5 sheet clips, 5 plastic plugs and 2 bed clamps.			
Patient Transfer System kit for I-400 (hospital)		b2011a	2 bed brackets, 1 transfer bar, 5 transfer straps and 2 transfer sheets.			
Patient Transfer System kit for		b2001c	2 bed posts, 1 transfer bar, 5 transfer straps, 1 mattress cover (non-slip), 5 plastic plugs, 5 sheet clips, 2 bed clamps.			

Exclusively at

I-400 (Home care)









TEXAS

BARTON CHAIR Features



Headrest Adjustable in height, depth, and angle for proper support and LIMITLESS number of positions.



Tilt-n-Space Designed for pressure relief and positioning. Moves position WITHOUT affecting body angle of patient.



Adjustable Leg & Footrest ADJUSTABLE to 90° angle and 6" for added length and comfort.



Anti-Bacterial Fabric Anti-bacterial antimicrobial, breathable material that can be EASILY CLEANED or DISINFECTED in seconds.



Patient Transfer System Accessories

Available for BC Models H-250, I-400, and I-700



FLORIDA GEORGIA

Refer a Patient (205) 631-8915 | Fax: (877) 863-9308